

living & co.

Foot & Back Massager



INSTRUCTION MANUAL

CL-218E

AFTER SALES SUPPORT

☎ (NZ) 0800 422 274

CONTENTS

Overview	2
Safety Instructions	2
Items Supplied and Product Overview	5
Operation procedure	7
Cleaning and Storage	8
Disposal	8
Technical Data.....	8
Troubleshooting.....	9
Warranty	10

Overview

- This device is designed to relieve muscle tension and provide relaxation. The device may be used to massage the feet or back muscles. It is not a replacement for medically required treatment.
- The device is for personal domestic use only and is not intended for commercial applications.
- The device should only be used as described in these instructions.
- The following are excluded from the warranty: all defects caused by improper handling, damage or unauthorised attempts at repair. This also applies to normal wear and tear.

Safety Instructions

- This appliance is not intended for use by persons (including children) with reduced physical, sensory or mental capabilities, or lack of experience and knowledge, unless they have been given supervision or instruction concerning use of the appliance by a person responsible for their safety.
- Children must be supervised to ensure that they do not play with the device.
- The device has a heating function. People who are sensitive to heat must be careful when using the heating function.
- Keep sharp and pointed objects away from the device to ensure that it is not damaged.
- The device must always be disconnected from the mains power when left unattended and before cleaning it.
- In the interest of your own safety, you should inspect the device for damage each time before you use it. It may only be used if it has no visible signs of damage!
- If the device, the mains cable, or the mains adapter of the device is damaged, it must be replaced by the manufacturer or customer service department or a similarly qualified person (e.g. specialist workshop) in order to avoid any hazards.



Health-Specific Instructions

- **Only use the device if you have no physical illnesses or ailments.**
- The following people may use the device **only after consultation with their doctor**:
 - people who wear a pacemaker or similar medical implants;
 - people with heart disease or who suffer from epilepsy;
 - people who have back complaints or spinal disorders;
 - people who suffer from osteomyelitis, bone inflammation (in particular if the spinal column is damaged) or from osteoporosis;
 - women who are pregnant or have just given birth;
 - people who suffer from circulation problems or thrombosis;
 - people who have undergone an operation for which the treatment is not yet complete;
 - people with physical abnormalities;
 - people who have a malignant tumour.
- If you are unsure whether the device may be used, consult a responsible doctor.
- Do not use the device near to or on injuries, swelling, skin irritations, infected skin, warts, fungal infections, abnormal changes in the skin or body parts with no feeling.
- The device must not be used for longer than 30 minutes continuously. Do not massage one group of muscles for longer than 15 minutes.
- If you experience an unpleasant sensation while you are using the device, stop performing the massage immediately!
- Do not fall asleep while you are using the device!



DANGER – Danger of Electric Shock from Water!

- Use and store only in closed rooms.
- Do not use the device in rooms with high humidity.
- Never immerse the device, mains cable or mains adapter in water or other liquids and ensure that they cannot fall into water or become wet.
- Should the device fall into water, switch off the power supply immediately. Do not attempt to pull the device out of the water while it is still connected to the mains power!
- Never touch the device, mains cable or mains adapter with wet hands when these components are connected to the mains power.



WARNING – Danger of Electric Shock from Damage and Inappropriate Use!

- The mains supply socket must be readily accessible so that the connection to the mains can be isolated quickly in the event of a fault or hazard.
- Only use suitable extension cables whose technical data is the same as that of the device.
- In the event of damage / faults (e.g. unusual noises, odour or emission of smoke), switch the device off immediately and pull the mains adapter out of the socket outlet.
- The device must be switched off when the mains adapter is pulled out or inserted into the plug socket.



WARNING – Danger of Injury!

- **WARNING – Danger of Suffocation!** Keep children and animals away from the packaging material.
- **Danger of strangulation!** Make sure that the mains cable is always kept out of the reach of small children and animals.
- The device should only be used and stored out of the reach of children and animals.
The mains cable must also be kept out of the reach of children and animals. Otherwise they might use the cable to pull the device down on themselves.
- **Danger of crushing!** Do not insert any fingers or other body parts into the device. Do not open the bottom zip during use.
- **Danger of tripping!** Make sure that the connected supply cable does not present a trip hazard.



WARNING – Danger of Fire!

- Do not use the device in rooms containing easily ignitable dust or poisonous and explosive fumes and substances in the air.
- Do not use the device under a blanket or similar item or in the vicinity of combustible material.
- In the event of a fire: Do not extinguish with water! Smother the flames with a fireproof blanket or a suitable fire extinguisher.

NOTE – Risk of Damage to Material and Property


- Do not lie or stand on the device.
- Protect the device from heat, naked flames, direct sunlight, extreme temperatures, moisture and impact.

Items Supplied and Product Overview



PLEASE NOTE!

- Never remove the rating plate and any warnings labels!
- Check the items supplied for completeness and possible transport damage. If you find any damage to the components, do not use them, contact our customer service department.
- Remove all packaging and protective coverings.

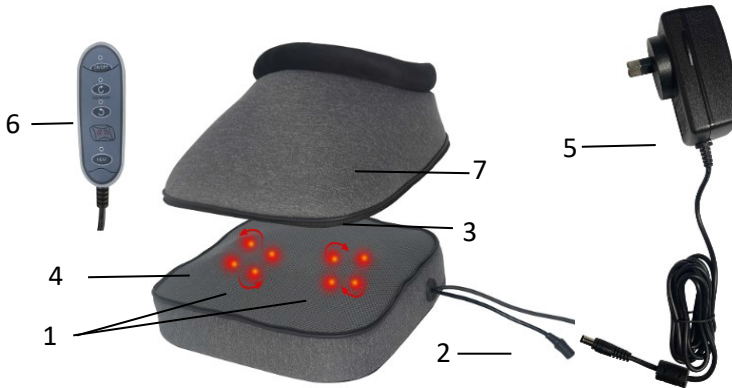


- 1 Massage nodes
- 2 Supply connecting cable / connector socket
- 3 Zip fastener (allowing attachment of upper foot enclosure)
- 4 Massage unit
- 5 Mains power adapter with output lead
- 6 Controller
- 7 Upper foot surround

NOTE:

Due to continued product improvement, images and illustrations in this manual may vary slightly from the product purchased. All images in this manual are for reference purposes only. Parts are not necessarily pictured to scale.

Controller

By pressing the buttons on the controller, the following functions can be chosen:

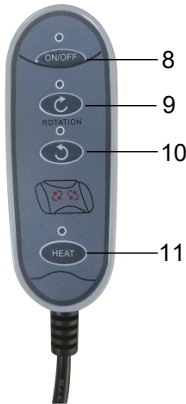

- 8** ON / OFF button
- 9** Clockwise rotation of massage nodes
- 10** Counter clockwise rotation of massage nodes
- 11** Heat / warming function on / off

NOTE: The activated function is indicated by the respective LED indicator lights above each control button.



PLEASE NOTE!

- No pressure should be exerted on the device as this will slow down or stop the massage and the device could be damaged. Do not stand or lie on the device.
- Do not use the device for longer than 15 minutes without stopping.
- For hygiene reasons, avoid direct contact of the skin with the device wherever possible. It should not come into contact with skin that has had cream or oil applied to it.
- A tingling sensation / itching during or after the massage is normal due to the increased blood flow circulation effect.

Operation procedure:

1. Insert the mains power adapter into an easily accessible mains socket outlet.
2. Connect the mains power adapter to the supply connecting cable / socket on the massager (2).



3. When used with the upper foot surround (7): Place your feet into the device and in direct contact with the massage nodes.
4. When used without the foot surround: Sit down on a chair with a backrest and place the massager so that it rests in direct contact with the desired area of your back.
5. Press the **ON / OFF** button (8) to switch on the device.
6. Press the desired button for Direction (9) or (10).
7. The massage will start to operate. To switch direction, simply press the corresponding button.
8. Press the Heat (11) button for the warming function to operate, this may take several minutes to generate an increase in temperature.
9. Press the Heat button at any time to switch off / on as desired.
10. After 15 minutes the device will automatically switch off. To switch off sooner, simply press the **ON / OFF** button once.
11. Disconnect the mains power adapter from the connecting cable / socket on the device.
12. Pull the mains power adapter out of the mains socket outlet and clean the device if necessary.

Cleaning and Storage



PLEASE NOTE!

- Do not use any caustic or abrasive cleaning agents for cleaning, these can cause permanent damage to the materials.
- Keep the device in a dry location at temperatures of between 0 and 40°C.
- Make sure that the device is switched off and the mains adaptor has been removed from the socket outlet.
- Wipe down the device with a damp cloth if necessary. Then wipe it down fully with a dry cloth.
- Keep the device in its respective packaging in a clean, dry place which is protected from dust and sunlight and is not accessible to children and animals.

Disposal



Dispose of the packaging material in an environmentally friendly manner and take it to a recycling center.



Dispose of the device in an environmentally friendly manner, it should not be added to general household waste. Dispose of it at a recycling center for used electrical and electronic devices.

Technical Data

Model number: CL-218E

Timer: 15 minutes

Power: max. 24 W

Mains Power Adaptor:

Model number: XH1200-2000WA

Input: 100–240 V ~ 50/60 Hz, 0.8 A

Output: 12V,DC, 2 A

Troubleshooting

If the device does not work properly, check whether you are able to rectify the problem yourself. If the problem cannot be solved with the steps that are listed, then you may need to contact an authorized appliance service centre.



Do not attempt to repair a defective electrical device yourself!

Problem	Possible cause / solution
The device does not work.	<ul style="list-style-type: none">• Has the ON/OFF button (8) been pressed correctly.• Is the mains power adapter inserted correctly into the socket outlet?• Is the socket outlet defective? Try another one.• Is the connecting cable connected correctly to Outlet cable of the adaptor?• The device has switched off automatically after operating for 15 minutes. Press the ON /OFF button again.• There is too much force or weight on the device (this may cause the massage motors to jam)
No heat is produced.	<ul style="list-style-type: none">• Has the heat function (11) been switched on?• Wait for 2 – 5 minutes for the heat to transfer
The device gets too warm.	<ul style="list-style-type: none">• Was it in operation several times in succession for 15 minutes? A maximum of 2 successive massage cycles should be applied. Remove the mains power adapter from the socket outlet and allow the device to cool down. Always disconnect the device from the mains supply when it is not in use.
The device makes noises.	<ul style="list-style-type: none">• Operating noises are normal as the massager incorporates mechanically driven parts. mechanics and not a defect.
Smoke emitted.	<ul style="list-style-type: none">• Pull the mains power adapter out of the socket outlet immediately and contact the customer service department.

living & co.



Warranty

We pride ourselves on producing a range of quality home appliances that are both packed with features, and completely reliable. We are so confident in our products, we pack them up with a 2 year warranty.

Now you too can relax knowing that you are covered.

Customer Helpline NZ: 0800 422 274

This product is covered by a 2 year warranty when

Accompanied by proof of purchase.

IM Version 1.0 – June,2021